
Massaggiatore contorno occhi  
Eye Wrinkle Eraser  
Masseur le contourne yeux  
Fallen-Rädler  
Massajeador arrugas contorno de ojos

## ITALIANO

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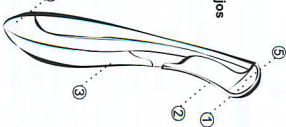## ENGLISH

### Product Feature

1. Intellectualized design, working automatically, no switch at all;
2. Iontophoresis and vibration massage function, anti-wrinkle effect of skin care essence or wrinkle care cream will be doubted at the help of this product. Introduce nutrition of skin care cream/essence into skin basal layer, refresh skin around eyes;
3. Introduce nutrition into skin basal layer, refresh skin around eyes; improve blood circulation and Lymph drainage, reduce eye puffiness. Eliminate eye pouch and remove wrinkle.
4. Mini & smart massage head, compatible with face T-Zone part (Lily);
5. Considerable protection cap as accessory to avoid massage head from dust and scratched.

## Cautions

1. This item use 1 YAAA battery. Please install correctly according to polarity mark in battery compartment. If battery counter-indialting, this product go into vibration/modification quality check mode) directly. It proved the product can be used normal.
2. Please don't use it in water.
3. Don't wash it in eyebath.
4. Clean the massage head with towel or tissue after use.
4. Please take out battery if not used for a long time.
5. Please keep out of reach of children.

## How to use

1. Slide out battery cover and insert 1 YAAA battery;
2. Apply desired skin care/eye care cream or essence on wrinkle part;
3. Hold this item keep silver plate touch, gently massage facial skin;
4. Please clean the massage head with towel or tissue after use, and put on protection cap.

## Product Structure

- 1 Massage Head
- 2 Cap

- 3 Battery Cover
- 4 Main Body
- 5 Sensor Plate

## TECHNICAL DATA

Voltage: 1\*YAAA alkaline battery  
Energy Consumption: 65mW  
Quiescent Current: <10µA  
Continuous Using Time: 7-9 hours

For any improvement reasons, Beper reserves the right to modify or improve the product without any notice.

## GUARANTEE CERTIFICATE

This appliance has been checked in the factory. From the date of original purchase a 24 month guarantee applies to material and production defects. The purchase receipt and the guarantee certificate must be submitted together in case of claim to guarantee.

CONTACT YOUR DISTRIBUTOR IN YOUR COUNTRY OR AFTER SALES DEPARTMENT BEPER.

E-MAIL: [assistenza@beper.com](mailto:assistenza@beper.com)

WHICH WILL FORWARD YOUR ENQUIRES TO YOUR DISTRIBUTOR.

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