LANAFORM®

Foot Tapping

Tapping massager

Massage by tapping

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LANAFORM FOOT TAPPING


Thank you for purchasing Foot Tapping by Lanaform®. This device is used to stimulate blood circulation, relax tense muscles, and relieve chronic fatigue.

Please read these instructions carefully before using this device, and please be certain to use it in a safe place.

Please read all instructions carefully before using this device.

NAMES AND FUNCTIONS OF THE PARTS

- 1) Indicator Light showing tapping strength
- 2) Indicator Light showing tapping speed
- 3) Indicator Light showing Mode of Usage
- 4) Mode-On/Off switch
- 5) Tapping Strength switch
- 6) Tapping Speed switch
- 7) Carrying handle
- 8) Tapping Cushion: The tapping motion can reduce tension in the sole of the foot
- 9) Power cord
- 10) Electric plug

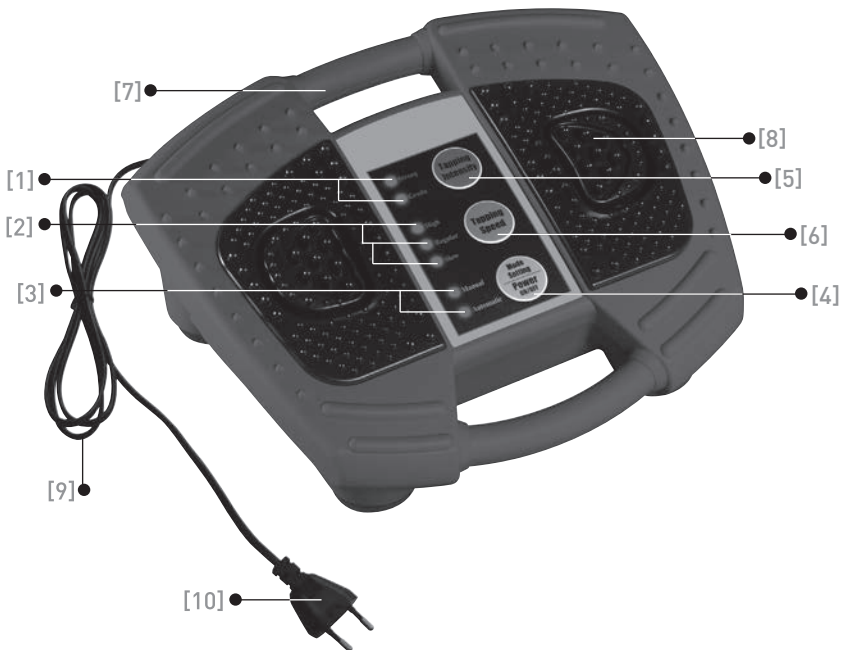INSTRUCTIONS FOR USE

Attention: Use of this massaging unit may interfere with television reception

- 1) Place the Foot Tapping on a stable surface and insert the plug into the electric outlet. Be sure that the Mode-On/Off switch functions correctly.
- 2) In a sitting position, place the tapping cushions over the tense or tired region, such as the sole of the foot or the calves. Concentrate on the soles of the feet to relieve muscle tension.
- 3) Press the Mode-On/Off switch. The Indicator Light will light up automatically. The timer is activated as soon as the power is switched on, and the massaging motion shuts off automatically after 15 minutes of usage.
- 4) When you press the Mode-On/Off button, the modes selected successively are: Power Off — Automatic Mode — Manual Mode — Power Off.

Automatic Mode: Modifies rhythm, tapping strength, and tapping speed automatically. Manual adjustments cannot be made to the tapping strength and speed.

Manual Mode: Tapping strength and speed can be chosen by the user.

- 5) When selecting Manual Mode by pressing the Mode-On/Off switch, the Tapping Speed and Tapping Strength switches will blink. Press the Tapping Strength and Tapping Speed buttons to adjust the strength and speed of the tapping according to your preference. To adjust the strength of tapping: pressing on the Tapping Strength button switches the indicator between High and Low. To adjust the speed of tapping: pressing on the Tapping Speed button switches the indicator between settings: High — Normal — Slow.

Attention:

- Do not block the tapping mechanism by pressing too hard on the cushion while massaging.
 - Do not use the cushion for longer than 15 minutes, the maximum period of usage. If you wish to continue massaging, please allow 30 minutes between each session.
 - Repeated and extended usage or too great a pressure on the tapping cushion could cause the internal temperature of the unit to rise and trigger the safety mechanism, stopping its motion. If this occurs, unplug the unit from the electric outlet and wait one hour before reusing the unit. If the unit remains plugged in, the tapping will restart immediately, posing possible danger.
- 6) The timer is activated when the power is turned on, and suspends the motion after a period of 15 minutes. If you stop the unit during a session, press the Mode-On/Off button and turn off the Mode indicator light. If you press the Mode-On/Off button, the unit will switch between modes: Off— Automatic Mode —Manual Mode — Off
 - 7) Unplug the power cord from the outlet after each use

MAINTENANCE AND STORAGE

Cleaning:

- Unplug the power cord from the outlet before handling. Not doing so could lead to electric shock or damage to the unit.
- Do not rinse the unit with water. The unit contains electric components and this could cause an accident or damage the unit.

Outside Casing of the Unit

- To clean the outside casing of the unit, dampen a soft cloth with a mild detergent diluted with water, wring it out, and wipe the outside of the unit.
- Do not use solvents such as paint thinner, benzene, or spray cleaning agents.
- Use cotton balls to remove dust on the surface of the cushion.

Storage:

Keep the device in its box in a dry and cool safe place and out of the reach of children. To avoid breakage, do NOT wrap the power cord around the unit. DO NOT HANG the unit from its cord.

Situations	What to Check	What to Do
There is no motion	Has the unit been unplugged?	Firmly insert the plug into the outlet.
	Is the electricity turned on?	Switch on the electricity.
	Is there damage to the power cord?	Contact your reseller.
Tapping motion stops during use	Has the 15 minute period elapsed since the unit was turned on?	The internal timer begins automatically and suspends movement after 15 minutes of usage. In this case, plug in the power cord again.
Suspension of tapping movement	Are your feet resting too heavily on the cushion?	If too much pressure is placed on the cushions or if the internal temperature has risen because of continuous use triggering the security mechanism which shuts off the unit, unplug the power cord from the outlet and allow the unit to cool down for one hour before reusing it.
The unit is damaged		Send it for repair.

IMPORTANT SAFETY INSTRUCTIONS

Danger:

to reduce the risk of electric shock:

- Always unplug the device immediately after use and before each cleaning.
- Do not try to pick up the device if it has fallen into water. Unplug it immediately.
- Do not use the device when taking a bath or shower.
- Do not place this device in any location where it could fall or be pulled into a bathtub or sink. Do not put it in water, or allow it to fall into water or into any other liquid.

Attention:

to reduce the risk of burns, fire, electric shock or physical abrasion:

- To avoid the risk of accidents or of any adverse reactions, anyone with any of the following symptoms or conditions should consult a physician before using this product:
 - Individuals under medical care or with physical abnormalities.
 - Individuals with a heart condition
 - Individual with a malignant tumour
 - Women who are in their menstrual cycle, or who are pregnant
 - Individuals who are under excessive stress or under duress
 - Individuals with back pain or an abnormality of the spinal column
 - Do not use this device in an area of the body that is currently receiving or has received medical treatment.
 - Individuals with a fever (above 38°C).
 - Individuals who are mentally impaired.
 - Individuals with fragile bones (for example, those who suffer from osteoporosis).

- The Foot Tapping should never be left unattended while plugged in. Unplug it from the electric outlet when not in use or before inserting or removing parts or attachments.
- Do not use the Foot Tapping with wet hands.
- Close supervision is necessary when this device is used with, or in the vicinity of children, or disabled or handicapped individuals.
- Please return the unit to the shop so that it may be returned to the LANAFORM service centre for inspection and repair.
- Do not use this device in any manner other than that prescribed, and only as described in this operation manual.
- Do not lie or sit on the Foot Tapping to obtain a massage as you could cause an accident, bodily harm, or damage the device.
- Do not use the Foot Tapping if the tapping cushion is not correctly installed. You could cause an accident, bodily harm, or damage the device.
- Do not use any device where the cord or plug has been damaged, or any device that does not operate correctly, has fallen or is damaged, or that has fallen into water. Please take it to a qualified technician or a service centre for inspection or repair.
- Do not disassemble the Foot Tapping or attempt to repair or modify it yourself. You may cause damage to the Foot Tapping and the unit will no longer be covered by the guarantee.
- Do not use the equipment near to a source of heat such as saucepans, heating, radiators, etc.
- Never use the Foot Tapping while sleeping. Do not fall asleep while using the massager.
- Never drop the Foot Tapping, and never insert objects into the openings.
- Do not use the Foot Tapping while spraying from an aerosol can, or while oxygen is being administered.
- Do not use the equipment near the eyes or any other highly sensitive areas.
- Do not use the equipment under bed covers or pillows. This could cause excessive heat leading to a fire, electric shock, or bodily injury.
- Do not move this device by pulling on the cord, and do not use the cord as a handle.
- Do not use the device outdoors.
- To unplug the device, switch it to the "Off" position, and then remove the plug from the outlet.

Attention

- Do not use this product for longer than 15 minutes. Avoid massaging the same area continuously. Elderly persons or those with a medical condition should begin with a short massage.
- Prolonged use could cause overheating, leading to a shortened life of the product. To prevent this from occurring, switch off the unit and allow it to cool before reuse.
- Never use this product directly on swollen or inflamed areas or in the area of a rash.
- This product is for non-professional use and is designed to provide a relaxing massage to tired muscles. Do not use this product as a replacement for medical care.
- Do not use this product before sleeping. Massage has a stimulating effect and could cause insomnia.
- This product should never be used by any individual suffering from physical problems that might interfere with their ability to operate the controls.
- In case of damage to the electrical cord, please take the unit for repair or servicing to a licensed technician in order to avoid any risk.
- This device has been designed only for use in the home.

TWO YEAR LIMITED GUARANTEE

LANAFORM® guarantees this product to have no material or manufacturing defect for a period of two years beginning from the day it was purchased, with the exception of the items below.

The guarantee of this LANAFORM® product does not cover damages arising from any abusive use or any incorrect use, accident, attachment of any unauthorised accessory, change made to the product or any other situation, whatever the nature, beyond the control of LANAFORM®.

LANAFORM® will not be held responsible for any incidental, consequential or special damage.

Any implied guarantees of product capacity are limited to a period of two years beginning from the initial purchase date.

Upon receipt, LANAFORM® will repair and replace, depending on the case, your FOOT TAPPING device, and will send it back to you at no charge. The guarantee may only be used through a LANAFORM® Service Centre. Any repair work on this product carried out by any person other than the LANAFORM® Service Centre will be deemed to annul the present

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Attention:

- Do not block the tapping mechanism by pressing too hard on the cushion while massaging.
 - Do not use the cushion for longer than 15 minutes, the maximum period of usage. If you wish to continue massaging, please allow 30 minutes between each session.
 - Repeated and extended usage or too great a pressure on the tapping cushion could cause the internal temperature of the unit to rise and trigger the safety mechanism, stopping its motion. If this occurs, unplug the unit from the electric outlet and wait one hour before reusing the unit. If the unit remains plugged in, the tapping will restart immediately, posing possible danger.
- 6) The timer is activated when the power is turned on, and suspends the motion after a period of 15 minutes. If you stop the unit during a session, press the Mode-On/Off button and turn off the Mode indicator light. If you press the Mode-On/Off button, the unit will switch between modes: Off — Automatic Mode — Manual Mode — Off
 - 7) Unplug the power cord from the outlet after each use

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Storage:

Keep the device in its box in a dry and cool safe place and out of the reach of children. To avoid breakage, do NOT wrap

the power cord around the unit. DO NOT HANG the unit from its cord.

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