



COLLAGENIUS®

V-FACELIFT AND NECK CONTOUR DEVICE

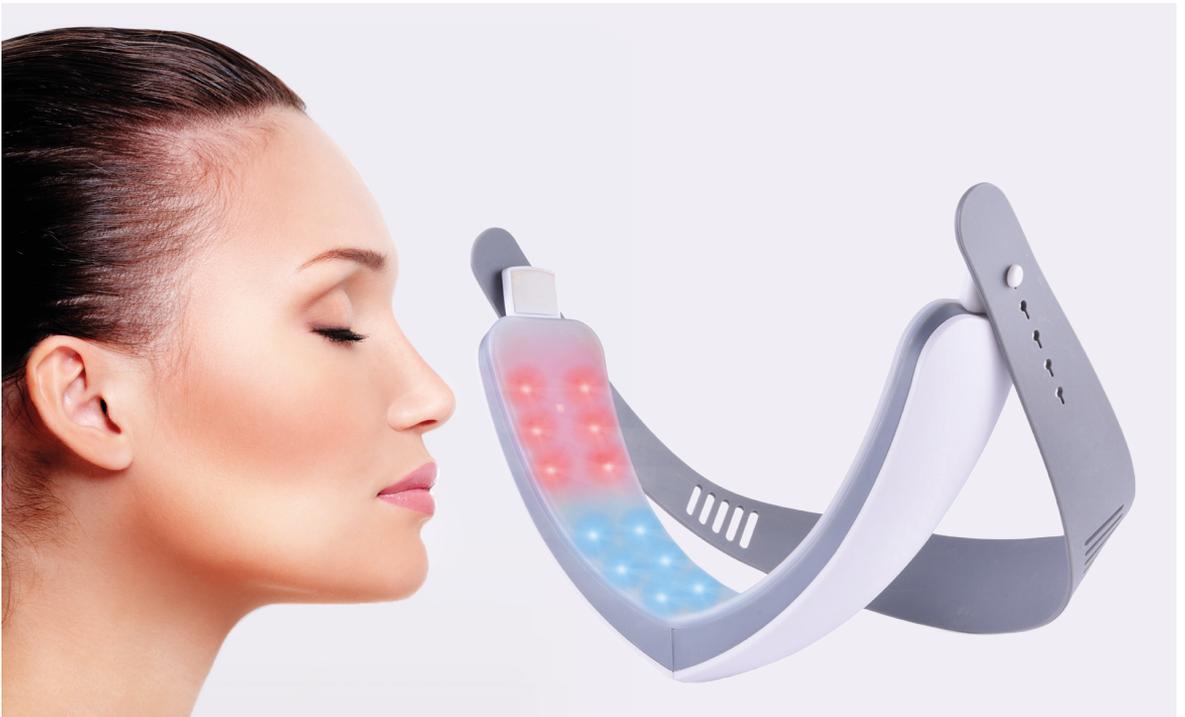

USER MANUAL

Please read this manual carefully,
keeping it for future reference.

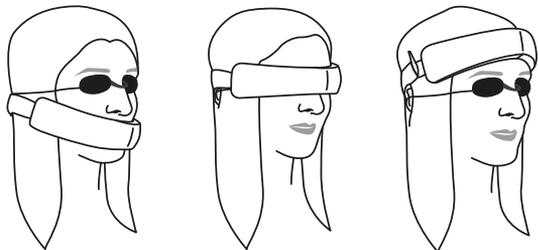

WELCOME

Thank you for choosing the Collagenius® V-Facelift and Neck Contour Device.

Your new Collagenius® V-Facelift and Neck Contour Device brings salon-worthy beauty into the comfort of your own home. Your device features three modes of colour light therapy alongside vibration therapy to:

- > Stimulate cell growth to clear your complexion
- > Improve circulation to revive tired, dull looking skin
- > Delay skin ageing and reduce fine lines and wrinkles
- > Tighten the skin to delay skin ageing
- > Improve muscle strength to contour the face
- > Tone the neck to reduce saggy jowls
- > Ease drooping around the nose and mouth



WHAT'S INCLUDED WITH YOUR COLLAGENIUS®

1. Collagenius® V-Facelift and Neck Contour Device
2. Adjustable band
3. Remote control
4. Protective goggles
5. USB charging cable
6. UK USB plug adapter
7. User manual



LED LIGHT AND VIBRATION THERAPY

Your Collagenius® V-Facelift and Neck Contour Device features 48 LED lights, offering superior LED light therapy alongside vibration therapy, to transform your skin and the underlying facial muscles.



RED LIGHT THERAPY to combat Ageing skin

Red light therapy helps improve the activity of skin cells, penetrating your skin at a wavelength of 625nm. Red light therapy promotes blood circulation, enhances and repairs skin elasticity, and contains anti-oxidant and anti-ageing properties. By increasing cell metabolism, your skin produces more collagen and fibrous tissue, enhancing a youthful complexion.



BLUE LIGHT THERAPY to combat Acne and Inflammation

Blue light therapy helps to calm your skin. Lightly penetrating at 470nm, it acts as an anti-inflammatory, regulating skin secretions and minimising acne damage. Blue light therapy works to clean and sterilise your skin's surface, balance natural oils, improve wrinkles and minimise pores.



INFRARED LIGHT THERAPY for a Glowing Complexion

Infrared light therapy penetrates your skin at 850nm. By increasing the temperature of your skin and expanding the micro vessels, the skin will be able to breathe more freely, and toxins are easily removed, leaving you with clear, acne-free skin.




VIBRATION THERAPY for a Contoured Appearance


Your Collagenius® V-Facelift and Neck Contour Device harnesses vibration technology to tone and tighten your face to give a more contoured appearance. The vibration therapy works at 6,500 revolutions per minute (rpm), helping to relax your facial muscles whilst also improving their strength and flexibility. Vibration therapy has been proven to improve the range of motion in the face, increasing blood flow under the skin.

USING YOUR COLLAGENIUS®


Use your Collagenius® 2-3 times a week to maintain a healthy, contoured appearance.

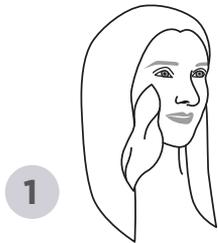


1 For best results, ensure the skin is free from makeup before use.

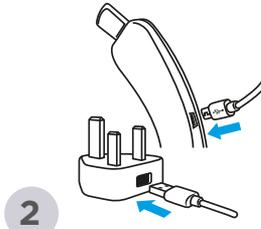


2 To charge your Collagenius®, plug the provided USB wire into the device and connect to a power supply. One full charge takes 150 minutes and lasts 300 uses.

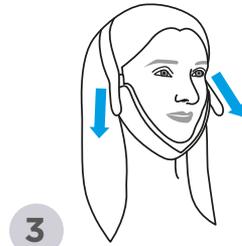


3 Attach the provided band to your Collagenius® and adjust to fit the shape of your face. See page 7 for ways you can wear your Collagenius®.

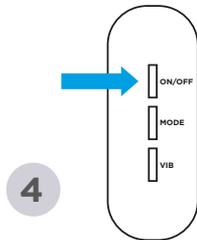


4 On the remote control, press the 'ON/OFF' button to turn on the power.

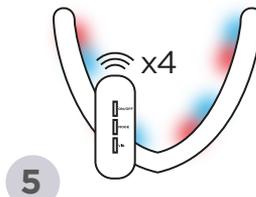

5 Your remote control will flash four times before automatically activating LED light and vibration therapy.

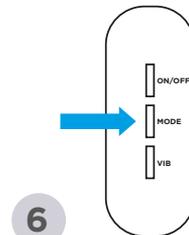

6 Press the 'MODE' button to select a programme mode (see page 6).

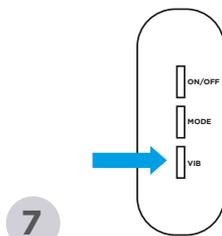

7 Press the 'VIB' button to change the vibration intensity. There are three levels of vibration intensity: low, medium and high.




8 After your ten minute session, your Collagenius® will automatically turn off. To manually turn off your Collagenius®, press the 'ON/OFF' button.

PROGRAMME MODES

Select the mode button to pick your favourite programme.


SKIN MAINTENANCE MODE

This mode alternates between vibration, red light therapy, infrared light therapy and blue light therapy, working to clear your complexion at the same time as contouring the muscles. This mode is recommended for optimum skincare maintenance, use 2-3 times a week.

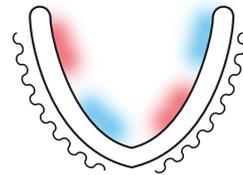


SKIN REPAIR MODE

In this mode red light therapy, infrared light therapy, blue light therapy, and vibration alternate. Vibration therapy works in harmony with the fast loop of light therapies for optimum skin repair. Ideal to clear acne and improve the appearance of pores. Use this mode as needed as one of your 2-3 sessions in a week.

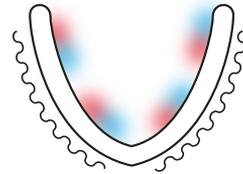


SKIN BRIGHTENING AND CALMING MODE

In this mode blue light therapy works exclusively to soothe the skin. This mode is recommended to relieve tired and dull looking skin. Use this mode as needed as one of your 2-3 sessions in a week.

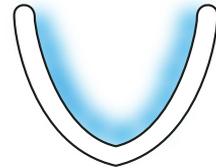

VIBRATION CONTROL

Use the VIB button to change the vibration to a comfortable level that tones your muscles for a youthful contoured appearance.

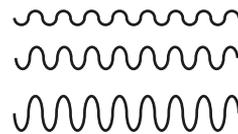

YOUR 10 MINUTE FACELIFT

Concentrate on specific areas of your face with the Collagenius® V-Facelift and Neck Contour Device.




PERFECT YOUR FOREHEAD

To clear imperfections, tighten the skin and reduce the appearance of fine lines and wrinkles, use the Skin Maintenance Mode on your Collagenius® across your forehead. To clear acne, use the Skin Relaxation Mode for a youthful glow.

SOLVE FINE LINES AND CROWS FEET AROUND YOUR EYES

Use the Skin Repair Mode on your Collagenius® over your eyes to boost collagen and reduce the appearance of fine lines.

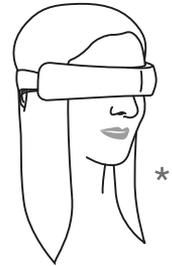

EASE DROOPING AROUND YOUR NOSE AND MOUTH

Use the Skin Maintenance Mode on your Collagenius® over your nose and lips to boost collagen production to tighten the skin and tone your muscles. Giving the appearance of a natural facelift.

PLUMP YOUR LIPS

Use the Skin Repair Mode on your Collagenius® over your lips to boost collagen production and reduce the appearance of fine lines around the lips. Giving you a naturally plumped look.



TIGHTEN LOOSE NECK SKIN AND REDUCE SAGGY JOWLS

Use the Skin Repair Mode on your Collagenius® to boost collagen production and tighten loose skin around the neck.

* We recommend using the protective goggles when using your Collagenius® in these ways.


CARING FOR YOUR COLLAGENIUS®

We recommend that you clean your Collagenius® V-Facelift and Neck Contour Device regularly. To correctly clean your device:

1. Turn it off and unplug from your USB charger.
2. Use a clean, dry cloth to gently wipe both the inside and the outer casing.

DO NOT use flammable liquids or harsh cleaners to clean this device.

After use, store in a cool, dry place. One full charge takes 150 minutes will last 50 hours. Do not leave charging for more than 12 hours. The remote control takes one 3V CR2025 battery.

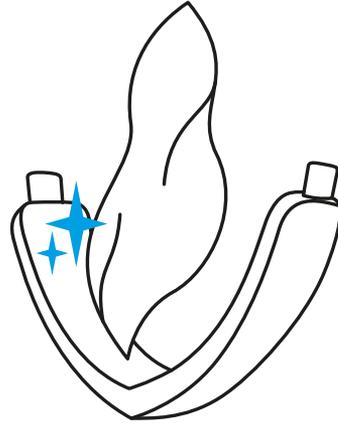

COLLAGENIUS®

Savantini House • Foster Street • Kingston Upon Hull • HU8 8BT • United Kingdom
Tel: 01482 888 785 • Email: info@savantini.com
www.stressnomore.co.uk

Safety Guidance and Contraindications

- DO NOT immerse the device or any of its parts in water.
- Use only as directed. If you experience any adverse reactions, discontinue use immediately.
- DO NOT use or store in temperatures above 35°C.
- DO NOT attempt to repair or modify this device. Doing so will nullify your guarantee.
- Keep away from children.

Do NOT use the device if:

- Your skin is severely inflamed.
- You are pregnant or breastfeeding.
- You have had facial plastic surgery.
- You have heart disease or complications.



Three Year Guarantee

A. This 3 year guarantee applies to the repair or replacement of your product, if found to be defective in material or workmanship. This guarantee does not apply to damage resulting from commercial, abusive, unreasonable use or supplemental damage. Defects that are the result of normal wear and tear will not be considered manufacturing defects under this guarantee. The manufacturer is not liable for incidental or consequential damages of any nature. This guarantee applies only to the original purchaser of this product from the original date of purchase. This guarantee does not cover damage to the adjustable band or USB charging cable. This does not affect your statutory rights.

B. Defective products should be returned to the place of purchase in accordance with policy.

C. This guarantee does not cover damage resulting from any unauthorised attempts to repair or from any use not in accordance with the instruction manual. This guarantee does not cover damages that attribute to self-dismantling and assembling, unauthorised operation not found in the instruction, random change of power source or voltage, artificial causes or natural disaster.