

IPL ProTM

Automatic Laser Hair Removal Machine

..... User Manual



IPL^{PRO}

AUTOMATIC FLASH & GLIDE

WELCOME

Thank you for choosing IPL PRO™ for your laser hair removal. IPL PRO™ offers salon-quality treatment in the comfort of your own home. Please read this user manual carefully before using your device.

YOUR IPL PRO™

Your IPL PRO™ Automatic Laser Hair Removal Machine is a home-use handheld, hair removal device that uses intense pulsed light (IPL) for long-term hair removal. IPL PRO™ can be used to remove unwanted hair from all over the body, including the bikini line, armpits, arms, legs and face.



The IPL technology used by IPL PRO™ provides long-term hair removal results at a fraction of the energy used by other light-based hair removal devices. Making it safe, as well as effective. IPL PRO™ protects your eyes and skin, designed so that a flash of IPL cannot be emitted unless the treatment window and skin sensors are in full contact with your skin.



HOW TO USE YOUR IPL PRO™

1. Shave, short-trim, epilate or wax the hairs you wish to remove. This ensures the IPL is absorbed deep into the hair follicles. See page 5 for guidance on using IPL to remove facial hair.
2. Clean your skin, making sure it is entirely dry and free from oils and creams.
3. Plug the power cable into a plug socket and press the power button to turn on your IPL PRO™. If this is your first time using your IPL PRO™, see page 4 for the skin test.
4. Choose your light intensity level by pressing the light intensity button.
5. Remove the IPL PRO™ device from the base unit and press the flash button once to remove it from standby mode.
6. Choose manual or automatic flash mode. See page 6 for guidance.
7. Place the IPL PRO™ device at a 90° angle against your skin surface. Ensuring the treatment window and skin sensors are fully against your skin.
8. If in **MANUAL mode**: Press and release the flash button to release a single flash of IPL. Glide your IPL PRO™ to the next area and repeat.

If in **AUTO mode**: Press the flash button once to remove it from standby mode, then press the treatment window against the skin and the device will release several flashes in a row. Slowly glide your IPL PRO™ across your treatment area.
9. After 30 seconds of non-use, your IPL PRO™ will enter standby mode. Simply press the flash button once to continue. To turn your IPL PRO™ off, hold the power button. It will automatically turn off after 10 minutes of non-use.
10. For best results clean your device after each hair removal session. Whilst unplugged, use a clean dry cloth to gently wipe the outer shell, skin sensors and treatment window.

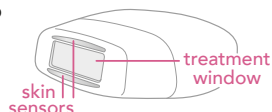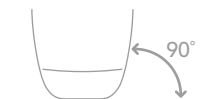


IMPORTANT - PLEASE READ

First Time Use - Skin Test and Choosing Your Light Intensity Level

The light intensity level setting determines the intensity of the IPL flash delivered to your skin. There are 3 levels – level 1 being the most suitable for fair skin, and level 3 for darker skin.

When you use your IPL PRO™ for the first time, or after sun exposure, perform a skin test on each area that you are going to treat to check your skin's reaction and to determine the correct light intensity level for each area.

1. Select light intensity level 1 and place the device at a 90° angle on the skin you wish to treat so that the treatment window and skin sensors are in full contact with your skin.
2. In 'MANUAL' mode, press the flash button to release a single flash of IPL.
3. Slide the device over the skin to the next area to be treated.
4. Increase the light intensity setting to level 2, and press the flash button again.
5. Repeat this for light intensity level 3 on a different area of skin.
6. Wait 30 minutes and check your reaction before proceeding with your treatment. Choose the highest setting that did not result in any skin reaction.

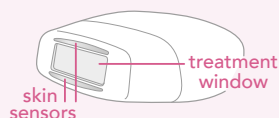

IPL hair removal may feel like a warm tingling sensation, but should never hurt. Do not increase the light intensity level if you experience any discomfort during or after use.

Choosing Manual or Automatic Flash Mode

Manual mode is ideal for areas that are small or curvy such as the knees and underarms. To activate manual mode, press the power button until the 'MANUAL' light is illuminated.

Automatic mode is ideal for large areas such as legs. To activate automatic mode, press the power button to switch to 'AUTO' mode.





Using Your IPL PRO™ on Your Face

If you wish to use your IPL PRO™ to remove facial hair, please always wear the eye protection goggles.

With the skin sensors on either side of the treatment window in full contact with your skin, you can use your IPL PRO™ on your:

- » Cheeks
- » Eyebrows
- » Moustache / upper lip
- » Sideburns
- » Jaw line

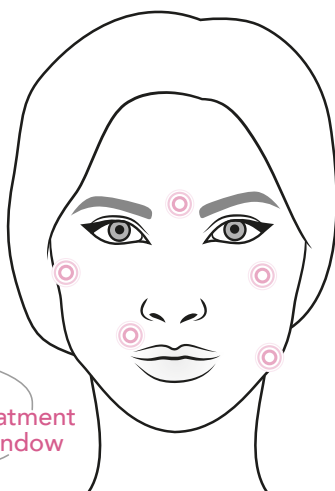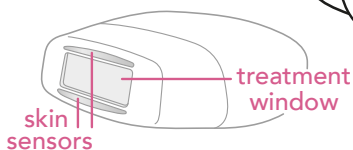

Please do not use IPL PRO™ for facial hair removal if your skin has had sun exposure or if you have undergone any professional facial treatment in the past 3 months.

SAFETY INFORMATION AND CONTRAINDICATIONS

Please ensure that you thoroughly read this information before using your IPL PRO™ Automatic Laser Hair Removal Machine.

1. Avoid using IPL on very dark skin pigmentations, including moles, large freckles, moles, warts, birthmarks or tattoos. This device may cause temporary pigmentation changes.
2. Do not use IPL after recent sun exposure or artificial tanning. Do not expose treatment area to sunlight for 1 hour after using the device.
3. If you naturally have white, gray or red hair IPL may not be effective.
4. Keep your IPL PRO™ away from water.
5. Do not use the device on the same area of skin more than once during each session and avoid overlapping flashes.
6. Stop using the device if treatment becomes uncomfortable. If your skin begins to blister or burn, stop using the device immediately.
7. Do not look directly at the light emitted from the device, wear protective goggles when using on the face.
8. IPL PRO™ offers permanent hair removal. Do not use the device on an area where you may want hair to grow in the future.
9. Never use flammable liquids such as alcohol (including perfumes, sanitizers, etc) or acetone to clean the skin before you use the device.
10. Do not cover the fans on the device or applicator.
11. Keep this device out of the reach of children.
12. If using a hair removal cream, wait 24 hours before using your IPL PRO™

DO NOT use IPL hair removal if:

1. You are pregnant or breastfeeding
2. You have sensitive skin, including eczema, psoriasis, lesions, open wounds or active infections. Wait for the area to heal before use.
3. You have a known light sensitivity, epilepsy, have a history of keloidal scar formation, are taking medication that makes the skin more sensitive to the light, or have a disease related to photosensitivity, such as porphyria, polymorphic light eruption, lupus etc.
4. You are currently or have recently been using alpha-hydroxy acids, beta-hydroxy acids, retin-a, accutane (isotretinoin), topical retinoids, azelaic acid or have taken steroids in the past 3 months.
5. You have an active implant, such as a pacemaker.
6. You have a history of skin cancer or potential skin malignancies or have received radiation therapy or chemotherapy in the past 3 months.

Please consult your GP or dermatologist if you are not sure if IPL is right for you.

Two Year Guarantee

A. This 2 year guarantee applies to the repair or replacement of your product, if found to be defective in material or workmanship. This guarantee does not apply to damage resulting from commercial, abusive, unreasonable use or supplemental damage. Defects that are the result of normal wear and tear will not be considered manufacturing defects under this guarantee. The manufacturer is not liable for incidental or consequential damages of any nature. This guarantee applies only to the original purchaser of this product from the original date of purchase. This does not affect your statutory rights.

B. Defective products should be returned to the place of purchase in accordance with policy.

C. This guarantee does not cover damage resulting from any unauthorised attempts to repair or from any use not in accordance with the instruction manual. This guarantee does not cover damages that attribute to self-dismantling and assembling, unauthorised operation not found in the instruction, random change of power source or voltage, artificial causes or natural disaster.



CONTACT US

ADDRESS

Savantini, Savantini House, Foster Street, Hull, HU8 8BT, UK

TELEPHONE

+44 (0)1482 496 931

(Available Mon - Fri, 8am - 6pm. Sat, 9am - 3:30pm)

EMAIL

sales@stressnomore.co.uk