


Eco Cosy® Back and Abdomen Heat Therapy Belt

With graphene and
far infrared technology

Please read this user manual before use.

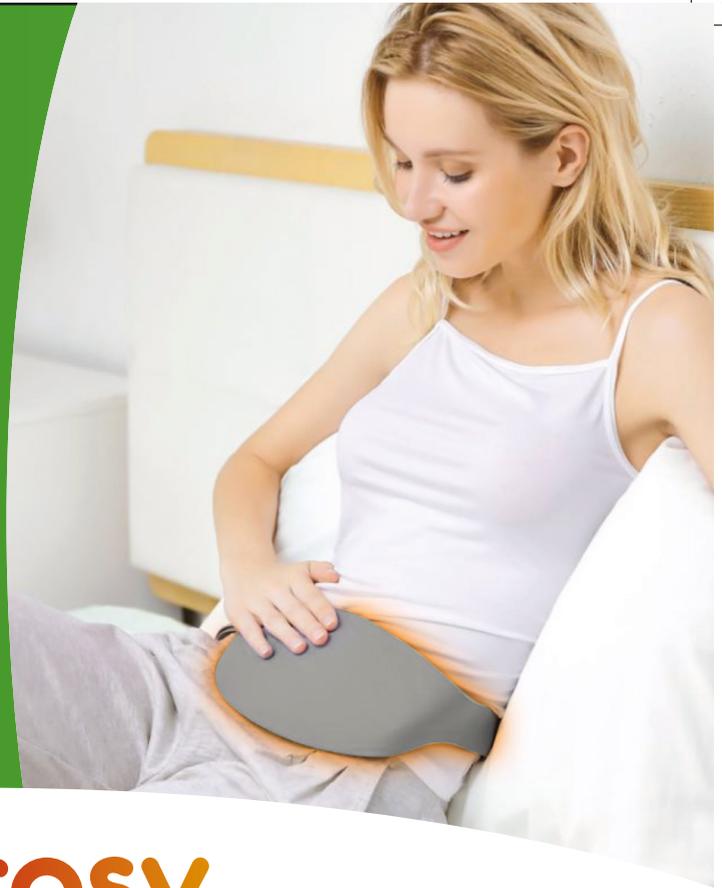

eco  **cosy**
warmth that won't cost the earth

Welcome

Thank you for choosing Eco Cosy®!

Eco Cosy® products create and insulate heat to help you reduce your heating bills and your impact on the environment. Whether it's for warmth and comfort, or much-needed pain relief, the Eco Cosy® Belt's intelligent temperature control technology uses graphene and far infrared to keep both you and the Earth comfortable.

With Eco Cosy® you won't have to deal with the guilt of turning the heating on. Costing just 70p per year* to use, compared to the average electricity bill of £699 per year, your Eco Cosy® Belt can save almost 1.5 tonnes** of CO₂ from entering the atmosphere each year.

*Based on 90 minutes of use a day with a standard 5000mha power bank (not included), and electricity costing 12p per kWh.

**Based on 3,200 kWh grid electricity used



Graphene + Far Infrared = Hi-Tech Eco Cosiness

Why Graphene?


Currently the world's lightest, thinnest and strongest new nano-material, graphene has a unique nano-structure that not only provides great physical flexibility, but also superior heat conductivity. Allowing your Eco Cosy® belt to heat up quickly.

Why Far Infrared?

Far infrared heat energy is of a similar wavelength to that which is radiated by the human body. Meaning it is absorbed deeply into the muscles, providing lasting pain relief and comfort. As well as improving blood circulation and reducing inflammation.



What's Included?

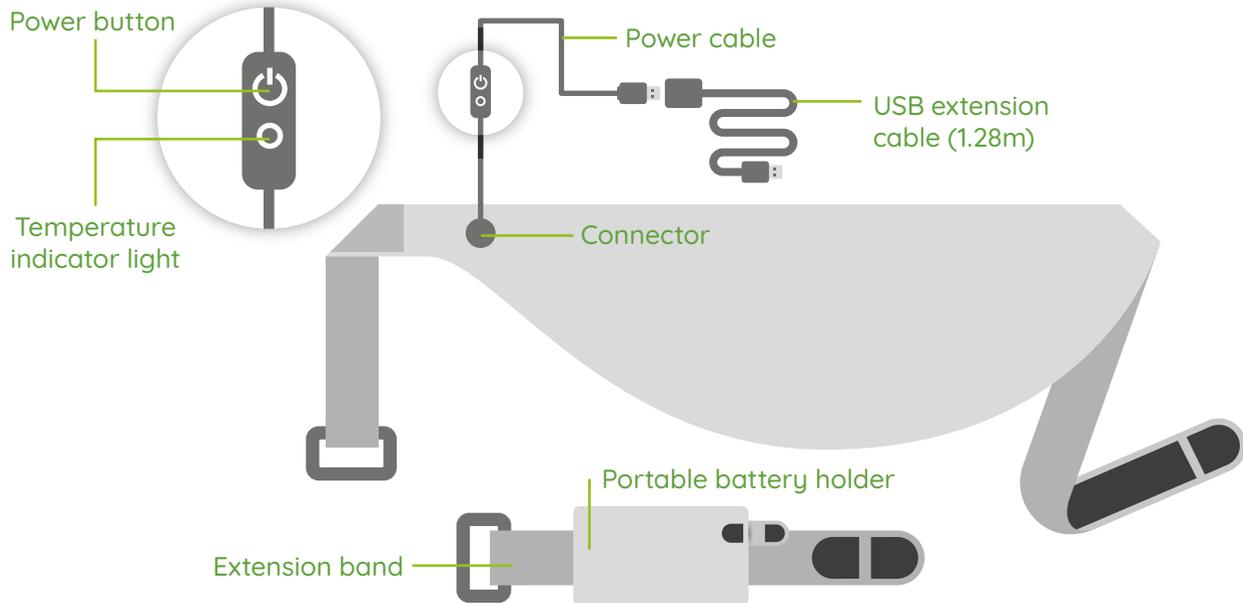

Instructions

1. Connect your Eco Cosy® Back and Abdomen Heat Therapy Belt to a power source using the USB extension cable. This could be a portable battery pack or a computer.
2. Hold the  button to turn your belt on.
3. Push the  button again to switch between temperature levels. The light below  indicates temperature level. Blue indicates the temperature is set to a low level. Green indicates the temperature is at a moderate level. A red light indicates the temperature is at a high level.
4. Wrap and secure the device around your waist, using the extension band if necessary.
5. When you are finished using your Eco Cosy® Belt, hold the  button to turn it off and disconnect it from the power.



Low
40°C



Moderate
50°C



High
55°C

Wearing Your Eco Cosy®

Your Eco Cosy® Back and Abdomen Far Infrared Heat Therapy Belt is slimline, breathable and can be worn discreetly under your clothes. You can even plug it into a portable battery pack for use on-the-go.


You can wear your Eco Cosy® Belt in a variety of positions.

Stomach - Stay warm in the cold weather or soothe an upset stomach by wearing your Eco Cosy® Belt around the stomach.

Back - Wear your Eco Cosy® Belt backwards to provide heated pain relief to your back.

Abdomen - Ease muscle aches and pains by wearing your Eco Cosy® Belt across your abdomen.

Lower abdomen - Soothe menstrual cramps across your lower abdomen with heat therapy from your Eco Cosy® Belt.

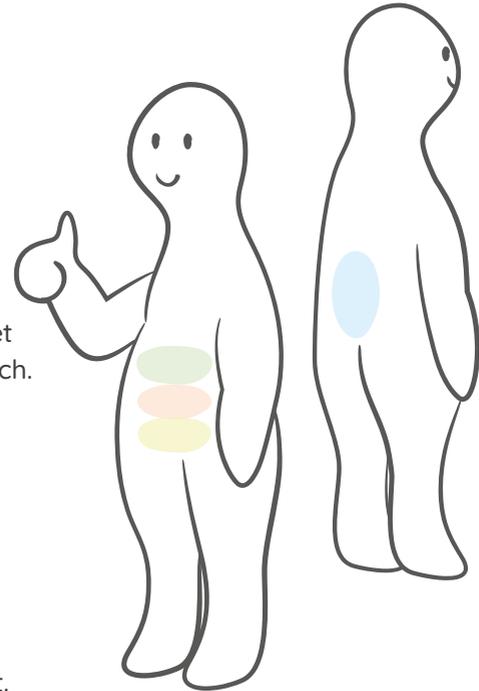

Precautions and Contraindications

- » Avoid using your device at the high temperature setting for prolonged periods.
- » When the device is not in use, ensure that it is turned off.
- » If wearing directly on the skin, lower the temperature setting to avoid discomfort and the risk of burns.
- » Do not use when sleeping.
- » If water infiltrates into the connector, stop use immediately. Wait for water to dry before reusing the device.
- » Do not disassemble or fold.
- » Do not allow children to use or play with this device.
- » Do not use if you are unable to use the device unaided, are hypersensitive, or suffer from gout.

Specifications

Power Input:	5V, 2A
Power Link:	USB
Power Rating:	5W
Service Environment:	-20°C to + 80°C Humidity ≤ 85RH%

Washing Instructions

1. Remove the power cable from the connector on your Eco Cosy® Belt.
2. Hand wash your Eco Cosy® Belt in cold water. Rubbing detergent onto any stains.
3. Lay flat on a towel to dry. Do not twist or fold.
4. Ensure the connector is completely dry before use.





Contact Us!

Email:

sales@stressnomore.co.uk

Phone:

+44 (0)1482 496 931

Our team of experts are

available at GMT:

Monday - Friday (8am - 6pm)

Saturday (9am - 3:30pm)

www.stressnomore.co.uk