

Folilaser

Folilaser Instructions Manual

Models:

Comb P1.0
Charger B1.0



"The industrial design of this product has been co-financed under the program "Innoempresa", as an innovation support for small and medium-large enterprises in the community of Madrid 2007 – 2013, by the Madrid Development Institute, Ministry of Economy and Treasury, Madrid Community, Ministry of Industry, Tourism and Trade, and by the European Regional Development Fund".

Index:

1. Scientific Principle of Folilaser
2. Description of Folilaser
3. System status
4. Use of Folilaser
5. Frequency of Use
6. Precautions
7. Guarantee

7. Guarantee:

Follaser guarantee covers any type of material or manufacturing flaws on the device for 2 years (24 months) starting the day of purchase and it has a one-year guarantee (12 months) -starting the day of purchase -on the accessories (such as the battery, power adapter, etc.) in normal application conditions. This guarantee is only issued to the original retail customer and only when the product was purchased at an authorized point of sale. If the product is defective within the guarantee period, please contact your authorized local provider.

Natural Logistics S.L. reserves the right to replace the defective product with a similar product in stock at the time of replacement.

This guarantee does not cover defective products due to:

- Accidents, inappropriate use, overuse or alteration
- Non-authorized technical service
- Use with non-authorized accessories
- Inadequate current or voltage plug in
- Any other uncontrolled situation

NATURAL LOGISTICS S.L. AND ITS AUTHORIZED PROVIDERS ARE NOT RESPONSIBLE FOR ACCIDENTAL, OCCASIONAL OR DIRECT DAMAGES THAT MAY RESULT FROM THE USE OF THE PRODUCT. ALL IMPLIED WARRANTIES, INCLUDING, AMONG OTHERS, IMPLIED WARRANTIES OF SUITABILITY AND MERCHANTABILITY, HAVE A 2-YEAR DURATION FROM PURCHASE DATE.

This warranty provides you with specific legal rights and there are also other rights depending on the country. Some States do not permit the exclusion of accidental or direct damages limitation or limitations on the implied warranty validity, so the previous limitations and exclusions may not be applicable to your case.

The Follaser laser comb is for domestic use and it has been professionally designed to prevent, correct and avoid problems related to hair loss.

Follaser offers a wide range of advantages that make it superior to other hair regeneration laser combs. This is mainly because it merges technology advantages with totally exclusive and registered designs boosting its efficiency while also providing the user with comfort due to its ergonomics –term previously unknown in this type of devices.

Please read this description and instructions carefully before using this product.

1. Scientific Principle of Follaser:

1.1 How will Follaser work on my hair?



Treatments with laser combs improve circulation, increasing the nutrients input to the cells and strengthening the tissues being treated. The application of these types of devices stops hair loss and stimulates hair growth and renewal together with its keratinization.

The light spectrum rays strike the skin increasing blood flow, oxygenation and nutrient provision to the tissues. Irradiation generates an increase in the tissues temperature and vasodilatation takes place enabling nutrients to better reach the tissues.

Blood stimulation to the scalp through the light wavelength appropriate to laser therapy of hair recovery enables the optimization of cell activity and allows a higher number of nutrients to reach the hair follicle and improve hair growth.

Due to the additional reconstructive properties of the skin, low power laser therapy can be used to heal the scalp after a hair transplant surgery.

Also, **Follaser** includes an innovative technology that cannot be found in any other comb in the market. This involves pressing a button to activate the vibration mode, VIBRATONIC, which makes the teeth vibrate and trigger the stimulation of the blood supply to the scalp, dramatically improving the results of laser therapy. Also, the upper button of the comb can control the vibration intensity according to each user preferences. There is no other hair recovery comb in the market offering such a state-of-the-art technology.

2. Description of the Product:

Follaser includes:

1. Laser Comb
2. Base Charger
3. Power adapter
4. 2 Interchangeable teeth of different hardness
5. Case

2.1 Laser Comb:

The comb includes:

Display



6. Precautions:

This device should not be used by people (including children) with limited physical, sensory or mental capacities or without experience or knowledge unless they have been supervised or trained on the use of this device by a responsible person in charge of their safety.

Only use this product for its original purpose. Do not use this product on other parts of your body. It should be used exclusively on the scalp. Children should be watched to make sure they do not play with the device.

This product is in conformation with EU Directive 2002/96/EC



The symbol of the "crossed trash can" in the product means that, at the end of its useful life, it should be disposed of separately from domestic waste. Therefore, it should be thrown away at a selective collection point of electric or electronic devices or it can also be returned to the provider at the time of purchase of a new similar product.

The user is responsible for delivering the device, at the end of its useful life, to the established collection points. The appropriate collection of the device allows the device to be recycled at the end of its useful life, its treatment and dismantling protecting the environment. This helps to avoid future negative effects on both the environment and general health and enables the recycle of the materials making up the product.

For more information on the available collection systems, contact your local collection offices or the providers that sold the product.



The device has an incorporated battery. Do not dispose of into regular waste after use. Dispose of in collection points indicated by local offices.

If the battery has been flat for a long period of time, an electrolyte leakage may occur damaging either the device or the user. Flat batteries should be removed from the device and disposed of at a collection point.

Do not open or brake the battery. The electrolyte leakage that could take place is corrosive and may cause burns or injuries to eyes or skin. It can even be poisonous if ingested.

Follilaser should be applied on dry scalp.


Four directions of application for the Laser Comb with Vibrator:

1° Stage: We start from the Frontal region all the way to the Occipital region of the head.

2° Stage: From the Occipital region to the Frontal region.

3° Stage: From Parietal to Parietal. From the left Parietal region to the right Parietal region.

4° Stage: From Parietal to Parietal. From the right Parietal region to the left Parietal region.

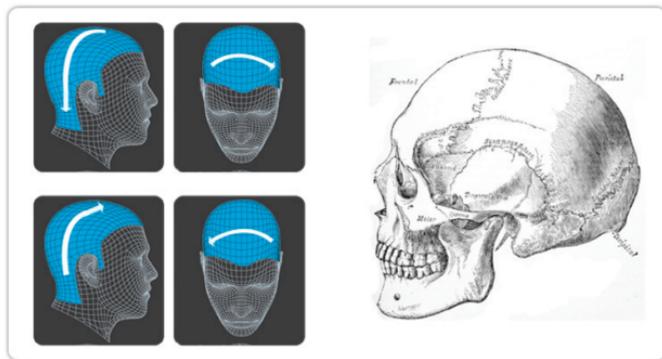

5. Frequency of Use:

The recommended frequency of application of Follilaser is 3 times a week but this may vary according to individual results. As a shocking treatment it is recommended to use Follilaser up to 4 times a week for the first 3 months so that the initial or recommended treatment -3 times per week- can follow next.

Likewise, the frequency can be decreased to 2 times a week -minimum- as a maintenance treatment as long as the recommended three-month treatment is over and the hair loss is not severe or just because the mere purpose is to keep a healthy hair.

The Display is constantly showing:

- **Remaining time of treatment. Countdown 15 minutes, shown in minutes and seconds.**

- **Work modes.**

Mode 0 (Default). No vibration. A zero (0) will appear.

Mode 1. Light vibration. Number one (1) will appear.

Mode 2. Strong vibration. Number two (2) will appear.

- **"Beep" activated: It will be active by default and it will show in the display. A speaker sign will be shown, if it disappears, it means the "BEEP" was deactivated.**

- **Charge Indicator: While charging, a battery will be shown at all times and the display will flash every 4 seconds.**

Button 1:

Press once and the treatment will begin. By default, the treatment will always last 15 minutes and will have a small "BEEP" activated that will sound every 4 seconds. To disconnect the "BEEP", only press Button 1 again and the "BEEP" sign will disappear from the display. To reactivate the "BEEP", Button 1 should be pressed again.

If, occasionally, we would like to turn off the device before the 15 minutes, Button 1 should be pressed for at least for 3 seconds. After 15 minutes, Follilaser will automatically turn off.

Button 2:

It will only work when Follilaser is on. Once on, Button 2 allows the user to choose the work mode.

- By default, if not chosen, the device will work in Mode 0 (No Vibration).

- If Button 2 is pressed once, the device will work in Mode 1 (Light Vibration).

- If Button 2 is pressed twice, the device will work in Mode 2 (Strong Vibration).

- If Button 2 is pressed three times, the device will return to its initial status, i.e., Mode 0 (No Vibration).

2.2 Base charger:



- Foilaser will be charged on this base charger. While charging, Foilaser will not turn on and it can only be turned on when the charger is not plugged in.

- When the charger is plugged in, it will show a red diode light on, which will not be lit otherwise.

2.3 Power adapter:

It is a power adapter of 9volts and 500mA.

Before plugging in the adapter, the pin should be plugged into the back part of the base charger. Always make sure the voltage matches the one indicated in the identification plate of the charger.

2.4 Interchangeable teeth:



- The teeth are interchangeable and there are 2 types of hardness to chose from. The one already placed is the softest but it can easily be switched by slightly pressing out with a finger on the upper joint. The upper part will come out first, and the lower part is held by 2 flanges that are easily removed when the teeth are moved upwards.

- To replace the teeth, the opposite process should take place.

3. System status:

3.1 Asleep:

In this work mode, LCD is off, in Power Saving Mode. Only Button  is the turn on button, which is the next step to ON status. If it is placed in the charger on this mode, it is automatically switched to CHARGING Mode.

3.2 On:

When Foilaser is activated, LCD will show the total time of treatment (15:00) that will go down per second.

To turn off Foilaser, keep the Turn On Button pressed for 3 seconds. Or otherwise it will do it automatically once the treatment time elapses.

The battery voltage is constantly checked so that if it is below the minimum value established (4.4V), the device automatically turns off. This way, if Foilaser is not charged and it's turned on, it will only remain on for 4 seconds showing low battery.

The smart activation of the vibration means that the system follows the same rhythm regardless of the battery charge so that the treatment is constant.

If Foilaser is ON and it is placed in the base charger, it is turned off automatically and it switched to CHARGING mode.

3.3 Charging:

This mode will show in the LCD by the battery symbol and the LCD backlight will turn on and off every 4 seconds.

Once the battery is charged to its maximum level, the LCD will turn off and stop the power input from the base charger.

If –at any work mode- the device is placed in the base charger and the battery is already charged (voltage over 5,5), the charge will not take place.

4. Personal Application of Foilaser:

Foilaser is used in a grid disposition on the scalp, acting on the 4 regions of the skull in 4 directions so that the laser diodes strike the scalp appropriately. The 15 minutes of the session will be divided equally among the 4 application directions -approximately 3 ½ minutes per direction at a medium-low speed.