# beurer

## FB 25


<b>D</b>	<b>Fuss-Sprudelbad</b> Gebrauchsanweisung.....	4
<b>GB</b>	<b>Foot Bubble Spa</b> Operating Instructions .....	7
<b>F</b>	<b>Bain de pieds bouillonnant</b> Mode d'emploi .....	11
<b>E</b>	<b>Baño de pies con masaje de burbujas</b> Modo de empleo.....	14
<b>I</b>	<b>Idromassaggio</b> Istruzioni per l'uso.....	18
<b>TR</b>	<b>Fişkırtmalı ayak banyosu</b> Kullanma Talimatı.....	21
<b>RUS</b>	<b>Гидромассажная ванна для ног</b> Инструкция по применению.....	24
<b>PL</b>	<b>Gorący masaż dla stóp</b> Instrukcja obsługi.....	28
<b>NL</b>	<b>Pedilúvio de hidromassagem</b> Gebruiksaanwijzing .....	31
<b>P</b>	<b>Pedilúvio de hidromassagem</b> Instruções de utilização .....	34
<b>GR</b>	<b>Ποδόλουτρο με αναβρασμό</b> Οδηγίες χρήσεως .....	37

Beurer GmbH • Söflinger Str. 218 • 89077 Ulm, Germany  
Tel.: +49 (0)731 / 39 89-144 • Fax: +49 (0)731 / 39 89-255  
www.beurer.com • Mail: kd@beurer.de

**CE IPX4**

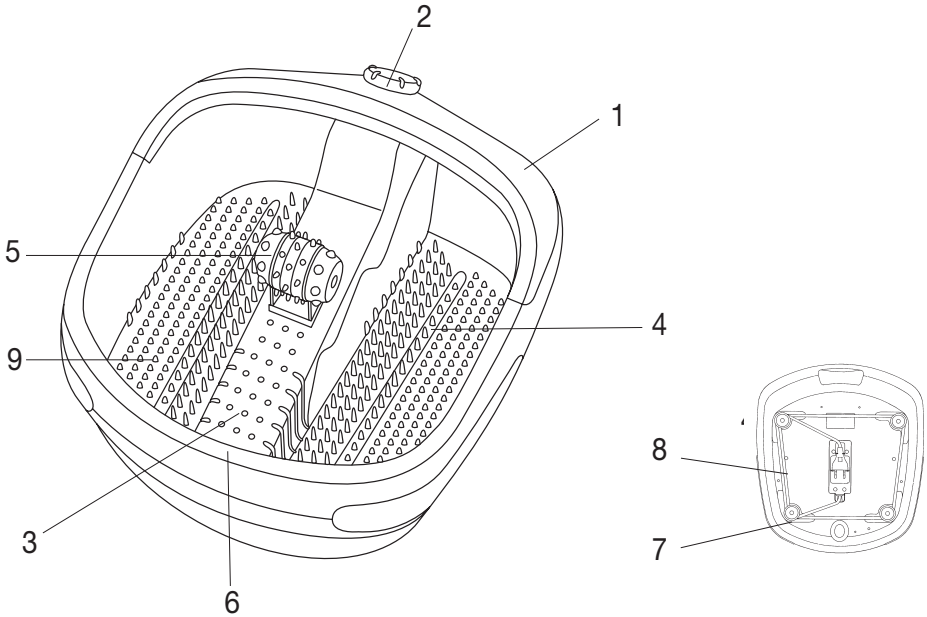**D Deutsch**

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## ENGLISH

### Contents

1. Beneficial for your feet.....	8	6. Important Information –	
2. Foot bubble bath FB 25.....	9	Keep for Later Use .....	10
3. Instructions for use.....	9	7. Technical specifications .....	10
4. Cleaning .....	9	8. Disposal.....	10
5. Exchangeable Massage Attachment for			
Massage of the Foot Reflex Zones.....	9		

### Dear Customer,

Thank you for choosing one of our products. Our name stands for high-quality, thoroughly tested products for the applications in the areas of heat, weight, blood pressure, body temperature, pulse, gentle therapy, massage and air.

Please read these instructions for use carefully and keep them for later use, be sure to make them accessible to other users and observe the information they contain.

Best regards,  
Your Beurer Team

### Items included in the package

- Foothbad
- These operating instructions

## **WARNING**

- The unit is only intended for domestic/private use, not for commercial use.
- This device may be used by children over the age of 8 and by people with reduced physical, sensory or mental skills or a lack of experience or knowledge, provided that they are supervised or have been instructed on how to use the device safely, and are fully aware of the consequent risks of use.
- Children must not play with the device.
- Cleaning and user maintenance must not be performed by children unless supervised.
- If the mains connection cable of this device is damaged, it must be replaced by the manufacturer, its customer service department or a similarly qualified person to avoid any risks.
- If water leaks from the appliance, the appliance should no longer be used.
- The device has a hot surface. Persons insensitive to heat must be careful when using the device.

## **Explanation of symbols**

The following symbols are used in the operating instructions and on the type plate:

**WARNING** Warning of risks of injury or health hazards.

**CAUTION** Safety information about possible damage to appliance/accessories.

**Note** Important information.

The device has double protective insulation and therefore complies with protection class 2.

Only use indoors

## **1. Beneficial for your feet**

Our feet play a weight-bearing role in our life – they carry our entire body weight. That is why your feet want to be looked after. There are various possible methods:

- **Magnetic field therapy:** the positive features of magnetic fields on the human body have been known for a long time. They are used in magnetic field therapy to increase the supply of energy to the cells.
- **Kneipp foot baths:** alternating hot and cold baths are the foundation of Sebastian Kneipp's health philosophy.
- **Foot reflex zone massage:** foot reflex zone massage is a very popular massage technique, in which the entire body is believed to be influenced by treating the feet.
- **Foot bubble bath:** the foot bubble bath is a simple and popular method for looking after your feet at home. Using a foot bubble bath, you yourself can actively do something good for your feet as often as you like.

## 2. Foot bubble bath FB 25

This foot bubble bath is extensively equipped:

- 16 magnets integrated in the foot bath
- Massage sole
- Detachable massage roller for foot reflex zone massage
- 3 functions: vibration massage, bubble massage and water heating

Relax yourself and your feet when using this foot bath. It can have the following beneficial effects:

- stimulates the circulation
- promotes cell regeneration
- revitalises and invigorates the feet
- relaxes tense muscles

### Note

warm foot baths with arnica, rosemary or sea salt extracts get the circulation going. Lukewarm foot baths with lavender or thyme have a relaxing effect and are a blessing for exhausted feet.

## 3. Instructions for use

Fill the foot bubble bath with warm or cold water up to the fill line (approx. 3 – 4 cm). Ensure that the rotary switch is set to "0" before the unit is connected to the mains.

Make yourself comfortable in a chair and place your feet on the foot rest and apply slight pressure. Never stand up in the foot bath. The rubber feet on the base of the unit prevent it from sliding.

Set the desired function using the rotary switch.

- Off (0)
- Heat and Vibration massage (W + M)
- Heat (W)
- Heat with Bubble and vibration massage (A + W + M)

A heating element is fitted below the foot surface so that hot water can be kept hot for some time. However, this heating element is not suitable for heating cold water.

The unit should not be used for longer than 40 minutes without interruption.

After your foot massage, turn the switch to "0" and remove the plug from the mains socket. The drain holes allow the water to be removed.

## 4. Cleaning

After use rinse the unit with tap water. The foot bath can be cleaned with a mild, non-foaming household cleaner, e.g. neutral vinegar. To dry the air lines, switch to the bubble setting without water for approx. 1 minute.

## 5. Exchangeable Massage Attachment for Massage of the Foot Reflex Zones

Additional comfort is provided by the local massage of the reflex zones. Massage of specific areas on the soles of the feet can exert an effect on organs of the body.

**It is recommended that you read the specialised literature and consult a doctor before the reflex zone massage to ensure that you obtain the desired effect.**

# 6. Important Information – Keep for Later Use

- **Read these instructions carefully and show them to other users.**
- Plug the unit only into the mains voltage specified on the rating plate.
- You must never connect the foot bath to the power supply when your feet are already in the water. If the device is damaged, this can lead to life-threatening electric shock.
- Ensure that the foot spa is undamaged before using it. If any fault or malfunction should occur during use, switch the unit off immediately and remove the plug from the mains socket. Never attempt to repair the unit yourself. Improper repairs can expose the user to considerable risks.
- The unit is intended only for the purpose described in the instructions for use.  
The manufacturer cannot be held responsible for damage caused by improper or careless use.  
The unit is not intended for commercial use, nor for use in medical clinics etc.
- The unit should not be used for longer than 40 minutes without interruption.
- Never operate the unit unattended.
- Never stand in the foot bath.
- Always switch the foot bath off and remove the plug from the mains socket before moving it or cleaning it.  
Always make sure your hands are dry.
- Never leave the device standing open over a prolonged period. Please use the relevant packaging for storage.
- Children are not aware of the dangers associated with electrical appliances. Make sure that the foot bath cannot be used by children without supervision.
- Never immerse the device in water.
- To fill the foot bath, it must not be placed in the shower, bath tub or wash basin.
- The device must not be connected to the power supply while being filled.
- Speak to your doctor: if you are experiencing pain in your leg or foot, e.g. phlebitis, varicose veins or diabetes, if you have a pacemaker, implants or other aids.
- The water must be drained only through the drain holes.
- Do not operate the foot massager without water.
- If water leaks from the appliance, the appliance should no longer be used.

## 7. Technical specifications

Mains connection 220-240 V ~ 50/60 Hz 130W

## 8. Disposal

Please dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment).

If you have any queries, please contact the appropriate local authorities.

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