


What's Included with Your Collagenius®

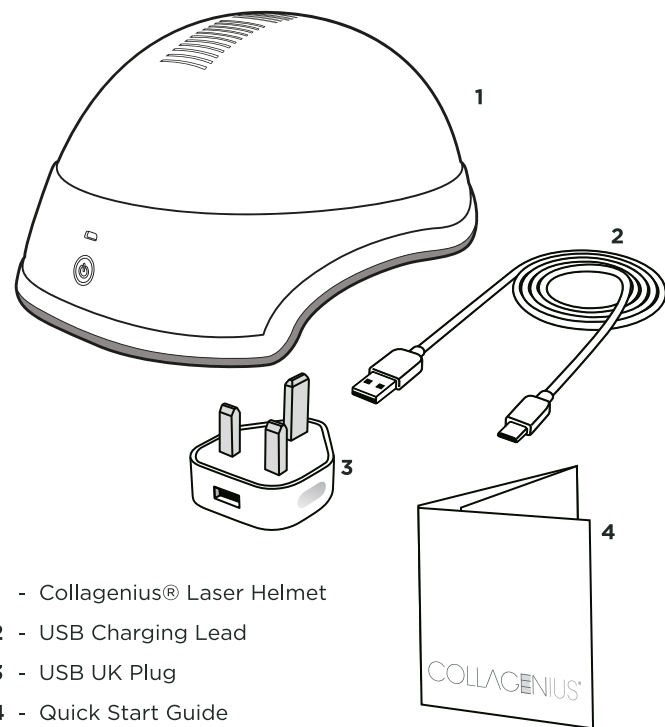

- 1 - Collagenius® Laser Helmet
- 2 - USB Charging Lead
- 3 - USB UK Plug
- 4 - Quick Start Guide

COLLAGENIUS®

Savantini House • Foster Street • Kingston upon Hull • HU8 8BT • UK

E-mail: info@savantini.com • www.stressnomore.co.uk



IFU12356_1221




COLLAGENIUS®

Hair Growth 144 Laser Helmet v2

Using Your Collagenius®

Step One

Before your first use your Collagenius® Hair Growth Laser 144 Helmet will need charging. A full charge will take around 4 and a half hours. The power button will illuminate red when charging and go off when complete.


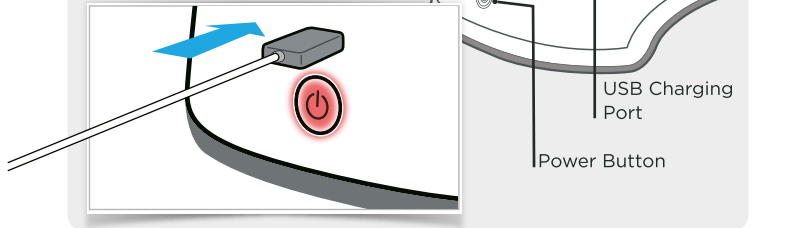


Step Two

Press the power button for 2 seconds, you will hear a 'beep' sound. The power button will turn blue. As a safety feature the laser lights will not illuminate until the helmet is placed on the head.


Step Three

Place the helmet on your head. We recommend wearing it with the power button at the back for comfort. The helmet laser lights will illuminate and start the 25 minute programme, to indicate it is working there is a 'beep' sound every 30 seconds.

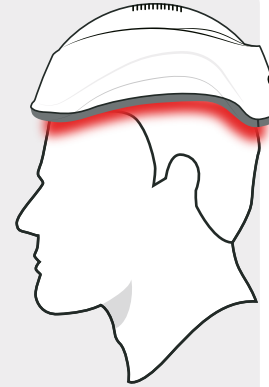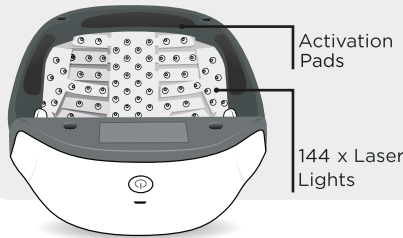


Step Four

If you wish to pause the programme simply remove it from your head. If you remove the helmet during the treatment time the laser lights will switch off for your safety. If you put the helmet back on within 6 minutes the programme will continue as normal. If the helmet is not put back on the head within 6 minutes, there will be a double 'beep' and will turn off automatically ending the programme.

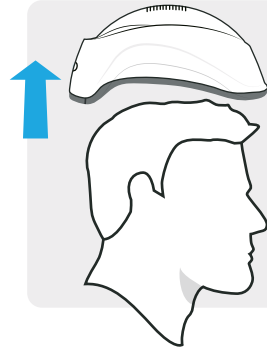

Step Five

At the end of the 25-minute treatment your helmet will switch off automatically. If you wish to switch off before then just hold the power button down for 2 seconds and you will hear a long beep. After use, store your laser helmet away safely in the box provided. It remains cool to the touch. One full charge should last between 4 - 6 treatments.

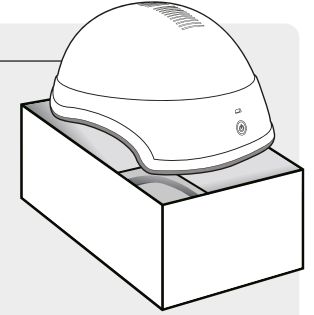

Frequently Asked Questions

How often should I use Collagenius®?

- Use every day for 25 minutes. You will notice significant results after 16 weeks, but you don't have to stop there, use it twice a week to continually improve growth and reduce hair loss.

Do I need to do any preparation before treatments?

- No, but your laser device should be used on a clean scalp so that the optimal laser energy can reach deep into the scalp for the best results.

Can I use the Collagenius® Laser Helmet with other hair growth products?

- Your Collagenius® device is completely safe to use with other hair growth products if you so wish to do so and won't reduce the effectiveness of the laser treatment.