



HYDRO 
GENER 

USER MANUAL

Please take the time to read this manual
before using your Hydrogener8™

THANK YOU FOR CHOOSING HYDROGENER8™

What is Hydrogener8™?

Hydrogener8™ produces hydrogen-rich water, providing unrivalled antioxidant care for your skin. Reversing environmental damage, reducing inflammation, and giving your skin a healthy glow.



YOUR ON-THE-GO SPA FACIAL

What is Ultrasonic Atomising Technology?

Your Hydrogener8™ harnesses the electrochemical production of ultrafine hydrogen rich droplets through a process known as ultrasonic atomising. These ultrafine droplets deeply nourish the skin. Quickly absorbing to refresh and rejuvenate.


What is Electrolysis?

Your Hydrogener8™ uses electrolysis to infuse water with added hydrogen. An electric current paired with vibration technology produces hydrogen rich water in a matter of seconds.

WHAT'S INCLUDED?

Your Hydrogener8™ Includes:

1. HYDROGENER8™ HYDROGEN WATER MIST SPRAYER
2. CHARGING CABLE
3. WATER REFILL BOTTLE
4. CARRY BAG

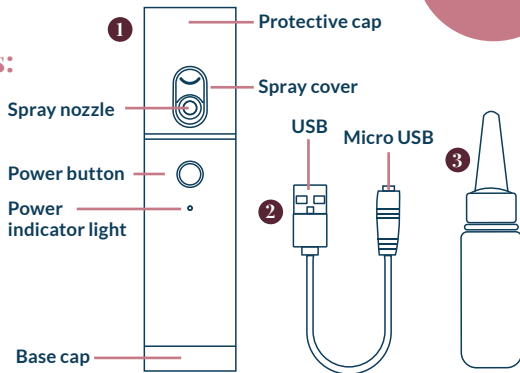


Charge by
USB or micro
USB

HOW TO USE YOUR HYDROGENER8™

Charging your Hydrogener8™

1. Remove the base cap and plug the charging cable into the bottom of the device (either USB or micro USB).
2. Connect the charging cable to a suitable power source. The power indicator light will flash while charging.
3. When fully charged, the power indicator light will stop flashing.
4. Unplug the charging cable and replace the base cap.
5. When the power indicator light flashes, recharge your Hydrogener8™

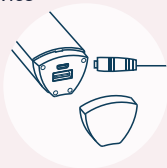

Filling your Hydrogener8™

1. Remove the protective cap by twisting it to the left and pulling upwards.
2. Open the rubber plug of the water tank and fill with mineral water using the provided water refill bottle.
3. Replace the rubber plug and protective cap.




One tank
of water can
last 10 uses

Using your Hydrogener8™

1. Open the spray cover to reveal the spray nozzle.
2. Hold the power button to turn on your Hydrogener8™.
3. A blue light will flash as your Hydrogener8™ begins electrolysis. After 10 seconds, the blue light will stay on to indicate that electrolysis has completed.
4. Hold your Hydrogener8™ approximately 10cm away from your face, press the power button once to release hydrogen-rich water mist for 30 seconds.
5. Move over the skin you wish to hydrate.
6. Repeat steps 3 - 5 as desired.
7. Hold the power button to turn off.
Close the spray cover.
8. Empty the tank when you are not using your Hydrogener8™ for a long period, and store in the provided carry bag.

Use under
or on top of
makeup for an
instant refresh

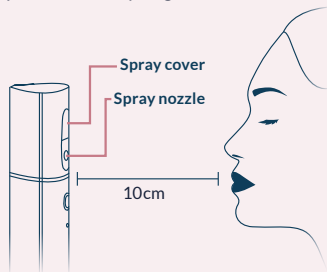

Safety Information and Contraindications

- For best results, use pure/mineral water
- Use the water refill bottle to fill up your water tank, avoid getting water elsewhere on the device.
- Do not use harsh chemicals to clean the device.
- Do not use on broken skin.
- Do not spray directly into eyes.
- Do not attempt to disassemble or modify the device.

Product Specifications:

MAX HYDROGEN WATER LEVELS: 1600 ppb

VOLUME OF WATER CONTAINER: 7ml

SIZE OF MIST: 0.3um

MIST FLOW: 0.6-0.8ml/30 seconds

USER TIME PER TREATMENT: 30 seconds

DIMENSION: 140 x 35mm

WORKING VOLTAGE: 100-240V 50/60HZ

INPUT VOLTAGE: DC5V/1A (Charging Cable)

OUTPUT VOLTAGE: DC5V/1A

CHARGING TIME: 3 hours

MAX POWER: 1.5W



www.stressnomore.co.uk

Savantini, Savantini House,
Foster Street, Hull, HU8 8BT

Contact our team of experts:

Mon-Fri: 8am - 6pm

Sat: 9am - 3:30pm

sales@stressnomore.co.uk

+44 (0)1482 496 931